

AE-248: Synthesis of [RuCl₂(CO)(p-cymene)]


Date: 2024-03-25

Tags: [RuCl₂(CO)(p-cymene)] NMR AE
reference analytics 1H reference
procedure

Status: Done

Created by: Alexander Eith

Reaction scheme/sample structure

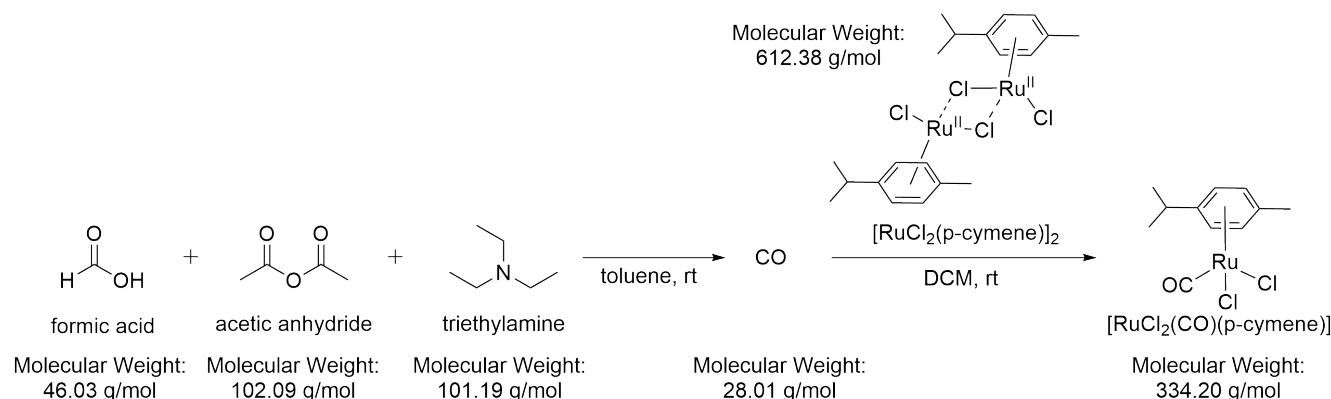

Literature/reference experiments

Literature	https://doi.org/10.1039/D1EE01053K
Reproduction	/
Related experiment	Experiment - AE-100: Synthesis of [RuCl ₂ (CO)(p-cymene)], Experiment - AE-221: Synthesis of [RuCl ₂ (CO)(p-cymene)]

Reagents

Name	CAS/Experiment number	Amount [mmol]	Equivalents	Mass _{theo} [mg]	Mass _{exp} [mg]	Molar mass [g/mol]	Volume _{theo} [mL]	Volume _{exp} [mL]	Density [g/mL]
[RuCl ₂ (p-cymene)] ₂	52462-29-0	1.30	0.5	800	794.24	612.38	/	/	/
Formic acid	Experiment - AE-247: Degassing of formic acid	2.60	1	119.6	/	46.03	0.098	0.098	1.22
Acetic anhydride Ac ₂ O	108-24-7	2.60	1	265.4	/	102.1	0.246	0.246	1.08
Triethylamine NEt ₃	Experiment - AEI-046: Degassing of NEt ₃ (Triethylamin)	5.20	2	526.24	/	101.2	0.72	0.72	0.73
DCM	Experiment - KRA-068: Degassing of DCM	/	/	/	/	/	16.0	16	/
Toluene	Experiment - AEI-016: Degassing of Toluene	/	/	/	/	/	5.2	5.0	/

Procedure/observations

All steps, unless mentioned otherwise, were carried out under argon atmosphere using standard Schlenk technique. (see Protocol [\[Protocol\] Schlenk Technique](#))

Date	Time	Procedure	Observation/Comments
25.03		The H-tube was prepared	
	9:20	The $[\text{RuCl}_2(\text{p-cymene})]_2$ was weighted at air in one of the two tubes (tube A). The flask was subsequently set under argon.	The $[\text{RuCl}_2(\text{p-cymene})]_2$ is a red glittering solid.
	9:45	DCM (16 mL) was added to tube A in argon counter stream.	A dark red solution was obtained.
	9:50	To the other tube (tube B) toluene (5.0 mL) was added in argon counter flow.	
	9:53	Fromic acid (98 μL) was added using a 100 μL Eppendorf Pipette in argon counterflow. The tip was flushed with argon for approx. 10 s prior to use. The flask was closed.	
	9:54	Acetic anhydride (246 μL , not dried, not degassed) was added using a 1000 μL Eppendorf Pipette in argon counterflow. The tip was flushed with argon for approx. 10 s prior to use. The flask was closed.	
		The reaction vessel was closed against the atmosphere and the Schlenk-line.	
	9:57 - 10:02	NEt_3 (0.72 mL) was added dropwise to tube B (approx. 5 min) while being stirred (520 rpm).	Some vapor formation during adding of first approx. 0.2 mL. After approx. 2 min the formation of gas bubbles was observed. Pictures: 20240326_095846.jpg 20240326_095916.jpg 20240326_100031.jpg 20240326_100056.jpg 20240326_100342.jpg 20240326_104117.jpg
	10:02 - 11:02	The reaction was stirred for 1:00 h at 520 rpm.	A dark red solution was obtained at the end. after reaction under CO.jpg
	11:02	The solution of tube A was transferred to the new 50 mL- schlenk tube (with a stirrer bar) using a syringe and a metal cannula. The flask was tilted so that non of the liquid in tube B could get into flask A, also the stop cock of flask B was kept closed	after transfer.jpg

	11:04 - 11:20	The solvent (DCM) was removed under reduced pressure using an external cooling trap.	after removing DCM.jpg
	11:21 - 15:45	The obtained solid was dried under reduced pressure for approx. 4:20 h.	A salmon red solid was obtained. dried solid.jpg
		The solid was stored under argon in the fumehood at rt.	
27.03	10:35	The solid was transferred into a vial and stored under argon at rt. During this transfer some material was transferred into a separately prepared 10 mL flask (approx. 5mg)	Significant amounts of product were lost due to the highly electrostatic behaviour
	12:20	To the flask containing to approx. 5 mg DCM-d2 (0.7 mL) was added	
		An NMR sample was prepared according to Protocol - Preparation of NMR Sample	NMR sample.jpg

Analysis

Date	Time	Sample name	Analysis method	Analytik device	Solvent	Raw Data	Processed Data	Comparative Data	Interpretation
27.03	14:16	AE-28-1	NMR 1 H	ZAF 300 MHz	DCM-d2	AE-28-1.zip	AE-28-1_10.nmrium	Experiment - AE-100: Synthesis of [RuCl2(CO)(p-cymene)]	analog to AE-100, some impurities are observed

Product characterization

	mass [mg]	purity [%]	mass _{pure} [mg]	amount [mmol]	Yield [%]	Description
AE-248-1	662.65	77 AE-248-NMR-purity.xlsx	510.24	1.53	58	salmon red solid

Future recommendations

Old procedure	Problem	Suggested new procedure
---------------	---------	-------------------------

Linked experiments

- [AEI-016: Degassing of Toluene](#)
- [AEI-032: Degassing of DCM](#)
- [AEI-046: Degassing of NEt3 \(Triethylamin\)](#)
- [AEI-057: Synthesis of \[RuCl2\(CO\)\(p-cymene\)\]](#)
- [AE-100: Synthesis of \[RuCl2\(CO\)\(p-cymene\)\]](#)
- [KRA-063: Degassing of formic acid \(80 %\)](#)
- [KRA-068: Degassing of DCM](#)
- [AE-221: Synthesis of \[RuCl2\(CO\)\(p-cymene\)\]](#)
- [AE-247: Degassing of formic acid](#)

Linked resources

- Protocol - [Preparation of NMR Sample](#)
- Protocol - [Schlenk Technique](#)

Attached files

- AEI-0013.cdxml
sha256: a5d83ef61c882538b0fa9abd8cd69e6c44ef56b9acd689ff66f798e4301229
- AE-248-NMR-purity.xlsx
sha256: 1367707f450e970502dae0382d71ab6d37683ef5a13685cbee1a2614b6f62299
- AE-28-1_10.nmrium
sha256: 0e23716ee8b3324373f447aba9d788c18317ebc06a714944108be1360eb7e0f3
- AE-28-1.zip
sha256: 73b84bf4637e90aaa5d200bdb697ed889aac2384a60824082d6510710d531964
- AEI-0012.cdxml
sha256: a5d83ef61c882538b0fa9abd8cd69e6c44ef56b9acd689ff66f798e4301229
- 20240326_100031.jpg
sha256: f9997283c2adec7dc88bf9d26e24c58d9f1c062a95329a4eefc080b12e806f52



NMR-sample.jpg

sha256: fbcac1765a017c63df59475f1d7c15bc11e568b593b5fecbd8ea78d09133b797



after-removing-DCM.jpg

sha256: c06833db009812b5384755dcaf88503ec302e3e385397a5f8e890d87afe9a6ae



20240326_095846.jpg

sha256: f7fc455d67d5e9f2cac1bee609eda71d88d4918525d082138ce1c1ac0d134463



20240326_095916.jpg

sha256: 8d14faffd3bf9f354f5cfcc4c6a4cc05d8426f2a4f872cf964a442063ee64955



after-transfer.jpg

sha256: 51d5ac4cf0d07c5e33b08026c77829e6a0e57984fbaf5511bc7f243923f7ffd9



20240326_100056.jpg

sha256: 538e55a0bb7ebbd087c43d12246ff138deb3018e8cfa787ec4682964c8b582a2



dried-solid.jpg

sha256: f5a72425ce771ef76ef8a5a54216945263d89aa2e2f65277e8a01af6279bd25f




20240326_104117.jpg

sha256: b5fa0e3bbcbd260a8a817f26a9a52d0b3f32429204812c1123c11d154bfaf960



after-reaction-under-CO.jpg

sha256: 1b77a01a67b666ec011e914128beadebd6ccf3f2abf7e40090578bcb8e5c7dc9

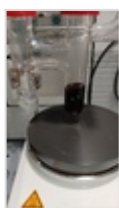

20240326_100342.jpg

sha256: fac549cd195674b2c0c78b95cad33e970cabcb02bf128a4c6b3aeb038c592c8



Comment

On 2024-10-25 09:00:40 Kristína Rabatinová wrote:

* missing the .png file and chemdraw file



Unique eLabID: 20240325-91b9b33765d7c8f75265b5c4c89c8cf7473697d2

Link: <https://elab.water-splitting.org/experiments.php?mode=view&id=914>