AE-261: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tertbutyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]


Date: 2024-04-18 **Tags:** [RuCl2(CO)(p-cymene)] Dichloro PNP [RuCl2(CO)(PNP)] P(tBu)N(py)P(tBu)


> NMR AE 1H 31P Status: Done

Created by: Alexander Eith

Reaction scheme/sample structure



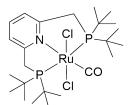

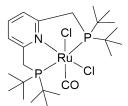

THF, reflux, o.n.

 $[RuCl_2(CO)(P^{tBu}N^{py}P^{tBu})]$

Chemical Formula: $C_{24}H_{43}Cl_2NOP_2Ru$

[RuCl₂(CO)(p-cymene)]

2,6-bis((di-tert-butylphosphaneyl)methyl)pyridine

Molecular Weight:

334.20 g/mol

Molecular Weight: 595,53152 Chemical Formula: C₂₃H₄₃NP₂ Molecular Weight: 395,55152

Literature/reference experiments

Literature	https://doi.org/10.1039/D1EE01053K
Reproduction	
Related experiments	Experiment - AE-157: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)] Experiment - AE-230: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)] Experiment - AE-232: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]

Reagents

Name	Amount [mmol]	Equivalen ts	Purity	Mass _{theo,p}	Mass _{theo} [mg]	Mass _{exp} [mg]	Molar mass [g/mol]	Volume [ml]
[RuCl ₂ (CO)(p-cymene)] (Experiment - AE-248: Synthesis of [RuCl2(CO)(p-cymene)]	1.53	1.00	0.77	505	657	1	334.20	/
PNP, 2,6-Bis(N,N-(tert-butyl)phosphinomethyl)pyridine (Experiment - AE-253: Synthesis of 2,6-Bis(P,P-di(tertbutyl)phosphinomethyl)pyridine	1.91	1.25	0.90	756	840	840	395.55	1
THF (Experiment - KRA-080: Degassing of THF	1	1	/	/	/	/	1	70
n-heptane (Experiment - AE-120: Degassing of n- Heptane	1	1	/	1	/	1	1	75

Procedure/observations

All steps, unless mentioned otherwise, were carried out under argon using standard Schlenk technique or in an argon filled glovebox. (see Protocol Protocol - Schlenk Technique)

Date	Time	Step	Observations		
18.04		A 150 mL Schlenk flask was prepared			
	14:40	The air cooler was dried by setting under vacuum	approx. 1:30 h		
	15:50	PNP(tBu) (840 mg) was weighted in the Glovebox K004 into the 150 mL Flask	white solid, PNPtBu.jpg		
	16:15	The whole product from AE-248: Synthesis of [RuCl2(CO)(p-cymene)] was used for the reaction	Ru.jpg		
	16:20	The flask which contained the [Ru] was washed with THF (2 * apprx. 10 mL) and the obtained solutions were transferred into the 150 mL flask	in-THF.jpg		
	- 16:30	To the reaction flask THF (50 mL) was added			
	16:35	The reaction mixture was stirred at rt and the flask was conected to the cooler in argon counterflow.	A red solution is obtained		
	16:40 - 17:10	The reaction micture was heated to reflux			
	17:10 -	The reaction mixture was heated to reflux for approx. 15:20 h at 350 rpm.	at-reflux.jpg after-reaction-at-reflux.jpg		
19.04	- 8:30	The reaction mixture was cooled down to r.t.	at-rt.jpg		
	10:15 - 10:35	The volume of the remaining reaction mixture was reduced under reduced pressure to approx. 35 mL	reduced-volume.jpg		
	11:15 - 11:30	To the solution n-Heptane (60 mL) was added under stirring.	after-adding-n-hept.jpg		
	13:10 - 13:45	The obtained mixture was filtered according to Filtration with frit technique and the residue was washed with n-Heptane (3 * 5 mL)	start-of-filtration.jpg		
	13:47 - 13:55	The obtained solid was dried under reduced pressure.	after first drying.jpg		

22.04	10:15 - 10:55	The obtained solid was dried under reduced pressure.	on fritjpg
	11:30	A NMR-sample was prepared from AE-261-1 according to [Protocol] Preparation of NMR Sample	NMR.jpg
26.04	10:30	The solid was transffered into a 10 mL Flask	AE-261-1

Analysis

Dat	e Time		Analysis method		Solvent	Raw Data	Processed Data	Comparative data	Interpretation
22.0	4 22:26	AE-261-1	1H, 31P (quant)	IAAC 400 MHz I	DCM-d2	AE-261-1_10.zip	AE-261-1_10.nmrium		Desired product at 65 ppm. Purity according to 31P of 94 %

Product characterization

	mass [mg]	purity [%]	mass _{pure} [mg]	amount [mmol]	Yield [%]	
AE-261-1	825.60	94	776.06	1.30	86	greenish solid, prod.jpg

Linked experiments

- JSC-KS-04: Degassing of n-Heptane
- AEI-047: Synthesis of 2,6-Bis(P,P-di(tertbutyl)phosphinomethyl)pyridine
- AEI-057: Synthesis of [RuCl2(CO)(p-cymene)]
- AEI-060: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]
- AEI-072: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]
- AE-100: Synthesis of [RuCl2(CO)(p-cymene)]
- AE-101: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]
- AE-102: Degassing of THF
- AE-108: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]
- AE-120: Degassing of n-Heptane
- AE-157: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]
- KRA-080: Degassing of THF
- AE-221: Synthesis of [RuCl2(CO)(p-cymene)]
- AE-230: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]
- AE-232: Synthesis of [RuCl2(CO)(PNP)] with 2,6-Bis((tert-butyl)phosphinomethyl)pyridine via [RuCl2(CO)(p-cymene)]
- AE-248: Synthesis of [RuCl2(CO)(p-cymene)]
- AE-253: Synthesis of 2,6-Bis(P,P-di(tertbutyl)phosphinomethyl)pyridine

Linked resources

Protocol - Preparation of NMR Sample

Protocol - Filtration with frit technique

Protocol - Schlenk Technique

Attached files

prod.jpg

sha256: d3b1e5ae30bd99a40fd678a2a35bbed02c35328809789b81629e873c95169215



NMR.jpg

sha256: a6e73b1966723ba79ebade8d7be02eb8a5626c144540e1f0edb495891eaeeb1c



on-fritipg

sha256: b5a3663697fc559e03dc6b8a4af078fbd245a4e16302b3d67b9e4fb95ed3fb18



AE-261-1 10.nmrium

sha256: d156fcb2ef6d59a9deeb7f85ebd86078c0f08e19084c909248092736dc4381e7

AE-261-1 10.zip

sha256: c41e58d463d8890db09ba18144754289a8d9adcd51306217ba39f35853a4488d

AE-1571.cdxml

sha256: 1f28659590507fe3cae942adc4ed5692c1008df067058eeaf4b8b10e69f9886e

Ru.jpg

sha256: c5ac21f0535891487ecfc7c8679c64fc8da3894f56a69874a10913b395f8b0ba



PNPtBu.jpg

sha256: e2bc3c9979b7947d6752351b701b88169b59cfe9c0fce6f2b6022588d716e6cf




after-first-drying.jpg

sha256: 264a6d96b3669ac657388c70adcd7f4701b46b7bf06faacfcae1f96085fbae38




at-rt.jpg

sha256: f643e79c8fce630c388d73896c476f911be2bd5e74db08bc483b7be4db11650f

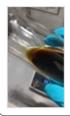

reduced-volume.jpg

sha256: 9fcbf809dd95387b62db89981401e5f8468a2af31db9df799084d0f8af656af9

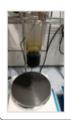

in-THF.jpg

sha256: 028679d0557932b31d97a0051f1b52ab18fa74e2606d2c8bb8531688d7b0eb8b



after-reaction-at-reflux.jpg

sha256: 0bc7960f70c01827c613c37c71ec06f30a6ee08914608a3c588b160231f956dd



after-adding-n-hept.jpg

sha256: cd9cc9879b3bf21d4199687938bf3715c9e7e0318a3f4ddbefc2b679cc283b6e



start-of-filtration.jpg

sha256: e88ae1390854cfe27956d613a4b928602cbadbffec393c464192cbabece4565d



at-reflux.jpg

sha256: 67812c9659cf049e3e40c00060cef8f800994ef4eed8b7d2579b6dec1b6eeba9



Comment

On 2025-02-19 15:08:13 Kristína Rabatinová wrote: Note: specify that whole ru storage schlenk was used



 $\label{thm:condition} \begin{tabular}{ll} Unique eLabID: 20240418-95059e9a21c50d0a979af71db3d34e36ff920a23 \\ Link: https://elab.water-splitting.org/experiments.php?mode=view&id=959 \\ \end{tabular}$