

Irradiation setup 1 (CEEC II, E004)

Date: 2025-09-29

Tags:

Lauda

Irradiation setup

EA

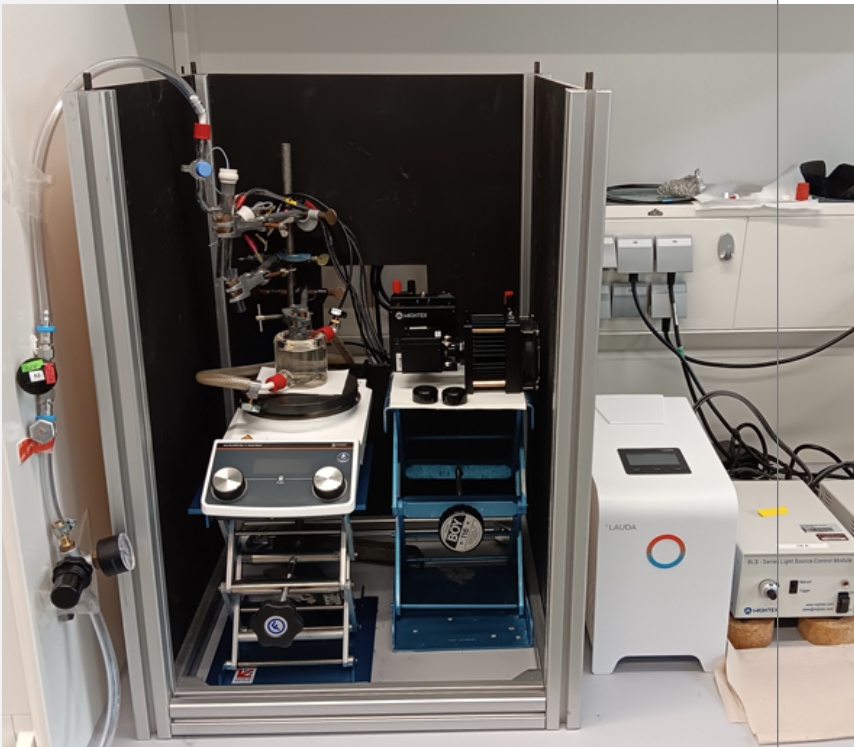
CEEC-II

LabE004

Category: Equipment

Created by: Ebrahim Abedini

Table

Setup name	Irradiation set-up 1	
Type of the setup	Equipment - Manual irradiation setup	
Information	<i>Information - Everything about Usage of Irradiation Set-Up</i>	
Location	CEEC II, lab E004 Lab E004 - CEEC II Left side of the lab Color: Yellow	
Name in Instruments Calendar	<i>Irradiation set-up 1</i>	
Power sources	Power Sources - BLS-18000-1 2 Power Sources - BLS-18000-1 3 Power Sources - BLS-13000-1E 2	
LEDs	Light Source - UHP LED 365 nm-2 Light Source - UHP LED 405 nm Light Source - UHP LED 6500 K (white)-2 Light Source - UHP LED 470 nm-2	
Beam combiner	Beam Combiner - Beam Combiner LCS-BC25-0409 -3	
Lauda	Equipment - Lauda - LOOP Thermoelectric thermostat, L100, 1	
FireSting / Sensor	Equipment - Firesting Fiber-Optic Oxygen Meter 2 Channel (Firesting 1) Equipment - OXSOLV-unprotected, QD7-576-089, CEEC II Lab E004 Equipment - OXSOLV-protected, QE7-580-089, CEEC II Lab E004 Equipment - Robust probe for liquid O2 measurement Equipment - Pt100 - Temperature Sensor	
Introduction	Alex, Nadja	
Responsible Person	Alex, Nadja	

Protocol on usage of the setup	Protocol - Irradiation setup Protocol - Catalytic tests in irradiation setup in flasks/vials
Additional Information	<i>Before the LED lamps can be turned on safety precautions have to be taken! (lab coat, blue nitril gloves, UV-glasses).</i> <i>Check the cleanness of the Lauda's solution before starting your photocatalytic experiments.</i> <i>Close the dark box from all 4 sides before starting the irradiation.</i> <i>The switch on the LED itself should not be used (should always be in the "on" ("I") position).</i> <i>After usage, check the Argon flow to be closed, LEDs disconnected from the power source, and the stirring plate to be turned off.</i> <i>Do not use the stirring plate's heating function while there are plastic tapes on it.</i> <i>Clean the degassing cannula after usage. (Replace with a new one in case of usage of noble metals).</i>
Status from	20/10/2025

Linked resources

Beam Combiner - [Beam Combiner LCS-BC25-0409 -3](#)

Equipment - [Firesting Fiber-Optic Oxygen Meter 2 Channel \(Firesting 1\)](#)

Equipment - [Robust probe for liquid O2 measurment](#)

Equipment - [Manual irradiation setup](#)

Equipment - [OXSOLV-protected, QE7-580-089, CEEC II Lab E004](#)

Equipment - [OXSOLV-unprotected, QD7-576-089, CEEC II Lab E004](#)

Equipment - [Lauda - LOOP Thermoelectric thermostat, L100, 1](#)

Equipment - [Pt100 - Temperature Sensor](#)

Laboratory - [Lab E004 - CEEC II](#)

Light Source - [UHP LED 365 nm-1](#)

Light Source - [UHP LED 405 nm](#)

Light Source - [UHP LED 365 nm-2](#)

Light Source - [UHP LED 6500 K \(white\)-2](#)

Light Source - [UHP LED 470 nm-2](#)

Power Sources - [BLS-18000-1 2](#)

Power Sources - [BLS-18000-1 3](#)

Power Sources - [BLS-13000-1E 2](#)

Protocol - [Catalytic tests in irradiation setup in flasks/vials](#)

Protocol - [Cleaning of double walled beaker](#)

Attached file

Irradiation-setup-1.png

sha256: a9ea29c9c81285ef836c11a372f8c8074e310c5de8a7d7b5f2571612a1995e42



Comment

On 2025-10-06 23:13:35 Ebrahim Abedini wrote:

* Pt100 added



Unique eLabID: 20250929-3fa8ac7f8618732062a286395cf106da242f3f12

Link: <https://elab.water-splitting.org/database.php?mode=view&id=279>