

Using the advanced irradiation setup V1.0

Date: 2024-11-18
Tags: Irradiation 3D print Irradiation setup advanced irrads setup
Category: Protocol
Created by: Alexander Eith

Goal

Irradiation using the Advanced irradiation setup V1.0

Prerequisites and preparation

Installed Advanced irradiation setup V1.0 ([Photoreactor - AE-367: Setup of advanced irradiation setup V1.0](#))

Holder and Holderdisk for desired reactor (e.g. 5 mL NS 14, NMR tube etc.)

Desired LED's (and their equipment) with power output characteristic

Aperture for desired LED

Be careful when handling the Holders, they can break or be deformed quite easily

Infos

The irradiation distance is 6.465 cm

Steps

Step number	Step description
1	Check the setup is complete and all parts you need are there as well
2	Prepare your sample
3	Install the desired aperture and the desired LED's (make sure, you know, which power setting you want to use) and secure the LED with a lab jack from falling down and moving
4	Take desired Holder and place it at the desired height (the flask should not touch the double walled beaker, but should also not be too high), the holder must be pushed as far as possible 20241118_112754.jpg
5	Check that the double walled beaker is filled with Milli-Q water, otherwise refill it
6	Take the desired holderdisk and place it onto the double walled beaker 20241118_112734.jpg
7	Place the reactor in the holder
8	Start your experiment
9	When the experiment is done remove the reactor

10	Finished, you can leave the setup as it is (make sure it is clean and that the LED's are standing on solid ground)
----	--


Linked experiment

Photoreactor - [AE-367: Setup of advanced irradiation setup V1.0](#)

Attached files

20241118_112734.jpg

sha256: 4affd79ff5da7c41c5b0b31f2525d202cbe1f98c4d2a8dd66d2f42e92c7b591f

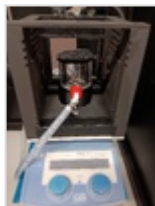

20241118_112754.jpg

sha256: 0c0e6a70db31df65c23a6ed316f85c1dd01774d03b81f6e872eafc7c130743d2



Unique eLabID: 20241118-85bc46733f4d211fd83221861ec1c9bba099f113

Link: <https://elab.water-splitting.org/database.php?mode=view&id=151>